
## Fusion System Anterior Matrix System

The all-new Fusion™ Anterior Matrix System from Garrison® gives you the materials and techniques to restore many types of anterior lesions with confidence.

For years, our customers have shown us how they use posterior metal matrices to restore anterior teeth. When placed vertically, what was the buccal-lingual curvature of the band creates ideal gingival-incisal curvature for an anterior tooth. We worked with these customers to design a metal matrix band that takes this technique to a whole new level of refinement.

Add to that the intelligent design of the new Fusion™ Anterior Wedges and you get a complete system. These radically curved wedges hug the bands to the tooth all the way from the facial aspect around to the lingual, helping to free up your hands. Composite placement and sculpting are simplified, allowing you to easily express your artistic nature restoring beautiful smiles for your patients.

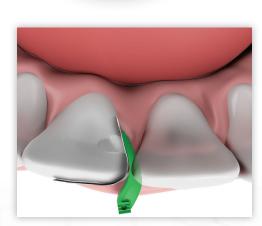

The Social Six never had it so good!



50 short bands

50 tall bands

25 extra small wedges

25 small wedges

25 medium wedges

25 large wedges







The **Fusion™ Anterior Matrix** contains a myriad of refinements to help you recreate ideal interproximal anatomy quickly and confidently.



## Fusion<sup>™</sup> Anterior Bands

		Standard			Bulk	
Refills	Item	Qty.	Price	Item	Qty.	Price
Short	AN100-M	50	\$78	AN100	100	\$137
Tall	AN200-M	50	\$78	AN200	100	\$137


- Only 0.0015" significantly thinner than strips to facilitate better interproximal contact.
- Firm stainless steel resists deformation both during placement and while compressing composite into the preparation.
- Easy to slide into the sulcus for deeper restorations.
- Ideal anatomic curvature in the gingival-incisal direction AND facial-lingual.
- Far superior to strips when restoring large or deep carious lesions.

The clever features of the **Fusion™ Anterior Wedge** make it a 'must have' for all anterior restorations.



## Fusion<sup>™</sup> Anterior Wedges

Description	Item	Qty.	Kit			
Fusion Anterior W	ANK4	400	\$172			
Fusion Anterior W	edge Kit			ANK4-M	200	\$90
		Standard			Bulk	
Refills	Item	Qty.	Price	Item	Qty.	Price
X-Small Yellow	ANYL-M	50	\$27	ANYL	100	\$49
Small Blue	ANBL-M	50	\$27	ANBL	100	\$49
Medium Orange	ANOR-M	50	\$27	ANOR	100	\$49
Large Green	ANGR-M	50	\$27	ANGR	100	\$49

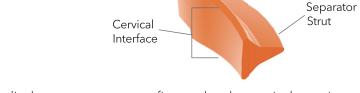

- Radical curvature ensures a firm seal at the cervical margin from facial to lingual.
- Helps maintain ideal anatomy by hugging the band to the tooth.
- Unique 'T' shape sits deeper in the interproximal space to minimize the dreaded 'black triangle.'
- Tapered separator strut allows you to adjust the amount of pressure exerted to control the interproximal contact.
- The cervical interface is half above and half below the sulcus for maximum adaptation with deeper restorations.

Be sure to check out the website at www.garrisondental.com for info on other innovative Garrison products, demo videos and our 6-month money-back guarantee.



150 DeWitt Lane Spring Lake, MI USA 49456 Offer expires 12/31/2021

To order: Call toll free 888.437.0032 or online at www.garrisondental.com





